Intelligent Face Mask

Quick Start Guide


Warning

This respirator helps protect against particulate contaminants but does not eliminate exposure to or the risk of contracting any disease or infection. Misuse may result in sickness or death. For proper usage refer to user instructions.

Disclaimer


JSB Tech Pte Ltd. does not assume any liability arising out of the application or use of any product; neither does it convey any license under its patent rights of others. In no event shall JSB Tech Pte Ltd. be liable for any direct, indirect, punitive, incidental, special consequential damages, to property or life, whatsoever arising out of or connected with the use or misuse of the product. Use the product by a user is the user's risk.

1. Rubber padding
2. Translucent face mask
3. Speaker module
4. Filter cap
5. Filters
6. Exhaust valve
7. Pump inlet cover

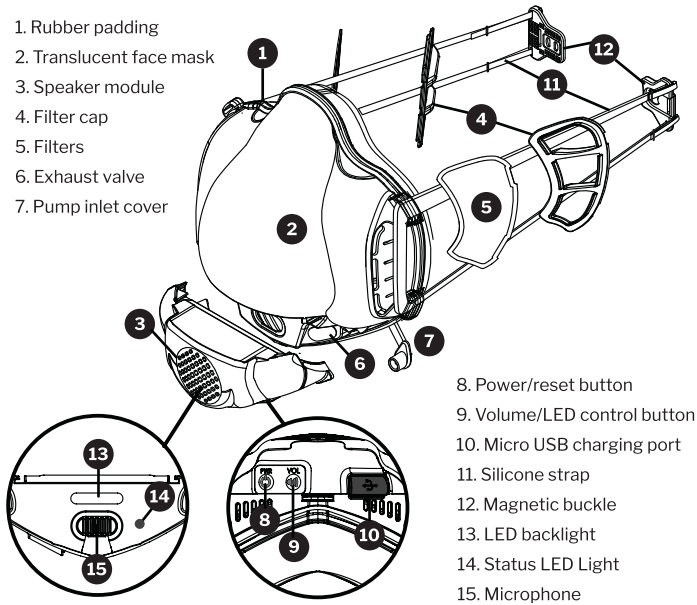How to setup & pair your smartphone:

Note: You will need to download "eGee-Pro Intelligent mask App" from Apple App store & Google Play.

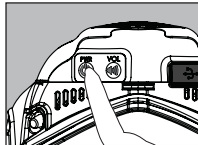Power on the mask by Press the "Power Button" located at the bottom speaker module.

The blue LED Backlight flashing, indicating is in **pairing mode**.

Tap on "Connect to Mask" button to begin the pairing.

JQC-6-09845-95-02

NOTICE

For better **COMFORT** and best **FITTING** on facial contour, please allow few days of wearing the mask to let the rubber padding and silicon straps soften.

Volume control

Users can adjust mask speaker volume up to 3 Levels, to change volume user press the volume button, each press will accompany by a 'short beeping' notification, with 3 beeps for max volume setting & further press will cycle back to level 1.

LED control

To turn LED light on/off just Press & Hold the 'Volume' button for 2 seconds.

Module Reset

To reset the module, just press & hold the 'Power' button for 4 seconds and release. The module will reboot automatically

NOTE:

User also can control the mask volume and LED colour through **eGee-Pro intelligent mask app**, available in Apple App store & Google Play. Scan QR below to download:

Using external USB to Charge the speaker module

Battery for speaker module can last over 3 hours* continuous usage after a full charge. Battery low notifications will appear in your mobile app or status LED will flashes in twice in green on the speaker module.


Use any standard external USB power source (e.g. a portable USB power bank) to power up & charge. When fully charged, the LED light will turns GREEN.

NOTE: * Battery life will vary depend on the volume setting & LED

How to remove / Install speaker module

To remove



Use both thumb press the microphone outwards from the face mask main body.

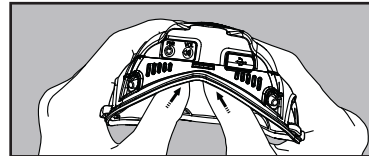To install

A) Align the module with the guide rail shown in the illustration and push toward the face mask until the user feels/heard a click.


B) Use both thumbs to hold the mask body in place and use index fingers and middle fingers to press the module towards the mask body untill hear/feel the click.

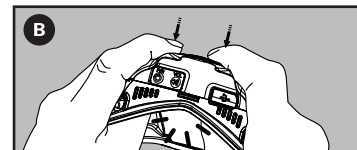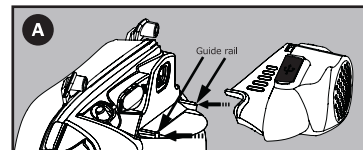How to prevent fogging

⚠ IMPORTANT:

Before washing & sanitizing the mask, remove all refill filters from the mask

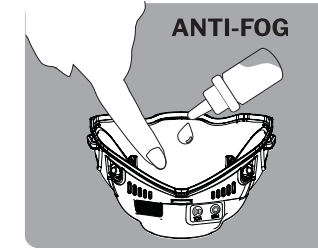- 1 Spray or wipe the soapy water on face mask and gently pat or rub dry using soft fabric.
- 2 Use a tiny drop of anti fog, or baby oil, or facial soap/gel, then spread evenly on the mask's inner surface using your index finger.

How to sanitize eGee-Pro Mask

CAUTION

The speaker module is not waterproof. When cleaning & sanitizing the mask, the speaker module must be detached from the mask.

Sanitizing Mask

Using soap/detergent sanitizing

Wash your mask in clean warm (around 30°C) soapy water, making sure your rinse all components thoroughly before allowing them to dry naturally

Sanitizing speaker module



UV light sanitizing

Place the module into UV light range for 15 Min, flip the mask to the front position and let it expose to UV light for another 15 minutes to sanitize.

Replacing filter

⚠ IMPORTANT:

- Properly wash & disinfect your hands before handling the facemask
- eGee-Pro antiviral refill filter service life is recommended for daily single usage and is limited by considerations on hygiene, damage to the mask, and increase in breathing resistance.

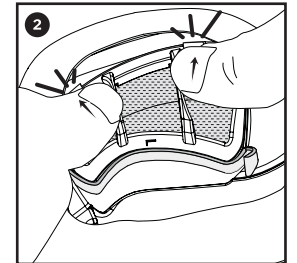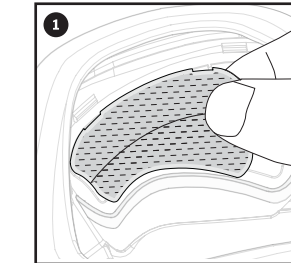1. Place the filter into filter gutter.
2. Close filter cap, push forward & press down untill you hear the clip sound.

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